
HX 2740



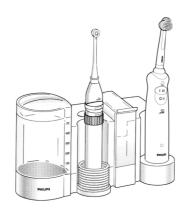

Philips Domestic Appliances and Personal Care

ServiceManual

PRODUCT INFORMATION

Type : HP 930
Description : Dental centre
Voltage : 220 - 240 V

: 200 V : 110 - 127 V

Power consumption : 20 Watts Voltage motor : 24 $\,\mathrm{V}$

SERVICE INFORMATION

Hand unit tooth brush (1)

Never open the hand unit of the tooth brush, because the unit is not water tight anymore.

NOTES

For standardisation reasons, the hand unit of the HX 2740 has to be replaced by the hand unit of the HX 2550.

The hand unit of the HX 2550 also has a charging light, that is green and blinks when the battery is being charged and that burns continuously when the battery has been fully charged.

When discarding the hand unit of the tooth brush, remove the built-in rechargeable battery.

The battery should be completely empty of charge.

Product meets the requirements regarding interference suppression on radio and T.V.

Published by Philips Domestic Appliances and Personal Care

Printed in the Netherlands

© Copyright reserved

Subject to modification





DISASSEMBLY ADVICE

PARTS LIST HX 2740

CO	N	N	E	C.	Т	N	G
	ΡI	F	C	F	п		

remove Wall holder I remove Screw A 2x remove Cover I 0 I

SWITCH 5

remove Wall holder I
remove Screw A 2x
remove Cover 101
remove Screw B and C
remove Cover 4
unsolder Wires switch

WATER JET 14

remove Wall holder I
remove Screw A 2x
remove Cover I 0 I
remove Connect.piece I I
lift Trafo/Mount.plate
remove Ring I 3/Tube

MOTOR 7

remove Wall holder I
remove Screw A 2x
remove Cover I 0 I
remove Cover 4
remove Connect.piece I I
unsolder Wires motor
lift Trafo/Mount.plate
remove Ring I 3 / Tube

TRAFO 8

remove Wall holder I
remove Screw A 2x
remove Cover I01
remove Cover 4
remove Cover 4
remove Connect.piece I I
unsolder Wires Motor
unsolder Wires switch
lift Trafo/Mount.plate

P.C.B. 13

remove Wall holder I
remove Screw A 2x
remove Cover I 02
remove Screw B and C
remove Cover 4
unsolder Wires


COIL 14

remove **Wall holder I** remove **Screw C 2 x** unsolder **Coil**

POS	DESCRIPTION	SERVICE CODE		
1	WALLHOLDER	4822 256 10189		
2	WATER CONTAINER	4822 441 81124		
3	CONTAINER BRUSH/SPOUTS	4822 418 10297		
4	COVER PRINTED	4822 442 01903		
6	MOUNTING PLATE	4822 464 30208		
7	PUMP UNIT	4822 360 10147		
8	TRAFO 220 - 240 V	4822 148 60303		
	TRAFO 110 - 127 V	4822 148 60304		
	TRAFO 220 V, KOREA	4822 146 10892		
	TRAFO 200 V, HONG KONG	4822 146 10893		
9	HOLDER WATER JET	4822 256 10188		
10	HOUSING PRINTED	4822 441 12374		
11	CONNECTING PIECE	4822 526 41098		
12	STUD	4822 462 40374		
13	P.C.B.	4822 214 70581		
14	COIL	4822 157 11413		
15	CLAMPING RING	4822 532 61262		
16	HAND UNIT WATER JET	4822 690 10227		
17	SET OF JET SPOUTS	4822 529 10439		
18	RECH. MOTOR UNIT NICd	4822 218 12093		
19	SET OF BRUSHES (2 IN BLISTER)	4822 479 20272		
Dir. for use	UK, D, F, NL, E, P, I, TR	4822 736 17096		
Dir. for use	Far East languages + UK	4822 736 17097		

IMPORTANT

- * AFTER THE PRODUCT HAS BEEN REPAIRED, IT SHOULD FUNCTION PROPERLY AND HAS TO MEET THE SAFETY REQUIREMENTS AND LEGAL REGULATIONS AS LAID DOWN AND ESTABLISHED AT THIS MOMENT.
- * PARTS, WITH NUMBERS FROM 100 UPWARDS, CANNOT BE ORDERED SEPARATELY.

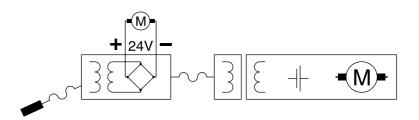

The Toothbrush charging is OFF when the Waterjet motor is ON

